



**virox**

# 70% IPA Disinfectant Liquid Spray

Effective bactericidal action against e.coli, pseudomonas aeruginosa and a diverse range of micro-organisms.



For more information, please call +65 - 6542 2100 or email to [enquiries@kleen-pak.com](mailto:enquiries@kleen-pak.com)

 Virox  [www.kleen-pak.com](http://www.kleen-pak.com)

### Overview

**Virox 70% IPA Disinfectant Liquid Spray** enables an instant and efficient disinfection to a variety of hard surfaces such as tables, countertops, wheel-chairs, surgery tables and x-ray equipment.

Ideal for all healthcare professionals who require an effective solution to their infection control needs.

Virox 70% IPA Disinfectant Liquid Spray is designed for hard surface disinfecting and cleaning and is durable for everyday use as well as non-smearing so they do not streak hard surfaces.

### Active Ingredients



Aqua, Isopropyl Alcohol 70% by Volume

### Antimicrobial Effectiveness Tested

Test	Method	Results (Percentage Kill)
<b>Antimicrobial Effectiveness Test</b>	<b>BS EN</b>	<b>Contact Time : 5 mins</b>
Escherichia coli	13727	>99.999%
Pseudomonas aeruginosa	13727/16615	>99.999%
Staphylococcus aureus	13727/16615	>99.999%
Enterococcus hirae	13727/16615	>99.999%
Canadida albicans	13624/16615	>99.99%
Asperbillus brasiliensis	13624	>99.99%
SARS-CoV-2	14476	>99.99%
Enveloped viruses	14476	>99.99%

### Specifications

Product Code	Volume	Bottle Spray Dimensions	EAN Barcode	Units / Carton	Carton Size & Weight	Pallet Weight	Pallet Height	Cartons / Pallet	Cartons / Layer	Layers / Pallet
VR4000	1L	D95mm x H290mm	8886389740001	6	L300 x W200 x H300 5.55kg	576	1650	100	20	5

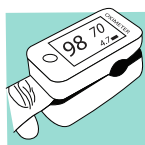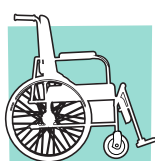


Suitable for disinfection of non-invasive medical devices in all areas and hard surfaces such as medical equipment surfaces, patient treatment areas, operating tables and similar work surfaces

### Features and Benefits



Kills 99.9% of SARS-CoV-2 & germs in 5 minutes

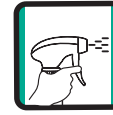

Designed for hard surface disinfecting and cleaning



This range carries CE Mark




Manufactured under ISO 9001, 22716 & 13485 quality system to ensure complete quality control




Formulated in Singapore

### Directions For Use

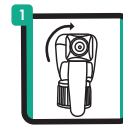

Turn nozzle to ON position

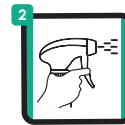

Keep the bottle upright, spray from 20-25cm away



Wipe over with clean damp cloth and allow to dry.